

HUWEL RT Mix

(Single Tube Reverse Transcription
Plus Real Time PCR Reagent)



Description

- HUWEL Single step RT Mix can be used for First-strand synthesis of cDNA from RNA molecules and Real time PCR in single tube.
- No need to prepare cDNA separately before setting up Real time PCR
- Just add RNA, primers and reagents to the tube and setup programme

Features

- Contains optimized mix A & B, to facilitate c-DNA synthesis real time PCR in single tube.
- Free from E.Coli DNA
- RNase and DNase free


Concentration

- 1 μ L /reaction of mixA and 12.5ul of mixB (for Final reaction volume of 25ul)

Applications

- First strand cDNA synthesis
- Single tube Reverse transcription Real Time PCR

Activity of HURT Mix

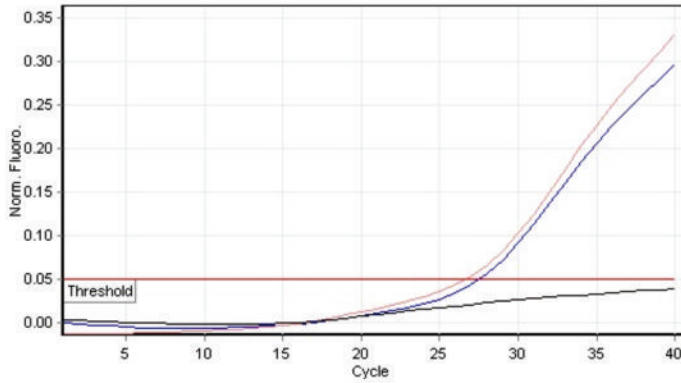

Serial	RT MIX	Sample Type	Ct.
1(Blue)	Commercial RT MIX	Standard	27.51
2(Red)	Huwel RT MIX	Standard	26.69
3(Black)	Huwel RT MIX	NTC	

Storage

- The recommended storage condition is -20°C .
- Avoid keeping at room temperature and repeated freeze thaws.

Ordering information

Cat. No	Reactions	Enzyme
HLRTM-25	25rxn	1 x 175 μl
HLRTM-50	50rxn	1 x 350 μl
HLRTM-100	100rxn	1 x 700 μl

Ref: Product and technical insert

OUR MOLECULAR DIAGNOSTIC KITS RANGE

- Quantiplus Adeno Virus Quantitation Kit
- Quantiplus BK Virus Quantitation Kit
- Quantiplus Epstein Bar Virus Quantitation Kit
- Quantiplus CMV Quantitation Kit
- Quantiplus Chikungunya detection Kit
- Quantiplus MTB/NTM detection Kit
- Quantiplus Dengue detection Kit
- Quantiplus HBV Quantitation Kit
- Quantiplus HSV 1&2 Detection Kit
- Quantiplus JAK-2 V617>F mutation detection Kit
- Quantiplus BCR-ABL Quantitation Kit
- Quantiplus Dengue genotyping
- Quantiplus JC virus Quantitation Kit
- Quantiplus JE virus detection Kit
- Quantiplus Cryptococcus detection Kit
- Quantiplus Malaria Detection kit
- Quantiplus Toxoplasma Detection kit
- Quantiplus HPV 16 & 18 detection kit

ALL KITS PACK SIZE: 25RXNS & 50RXNS



HUWEL LIFESCIENCES PVT. LTD.

Thriving with science

First Floor, H.No. 3-14/2, Survey No. 312, Narsingi Village, Gandipet Mdl, R.R. District, Hyderabad - 500 075, Telangana. • Mobile: +91 9052226226

E-mail : info@huwellifesciences.in • Website : www.huwellifesciences.in